


WILAméd

Equipment for Professionals

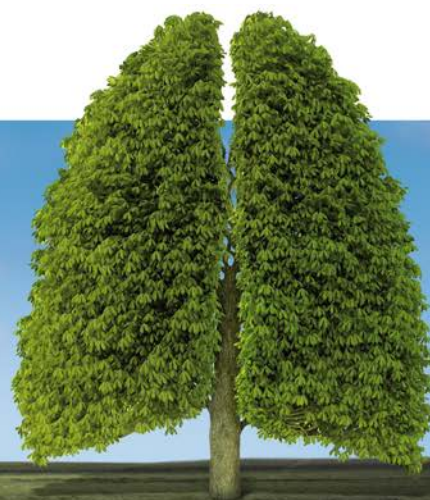

WILAflow Elite Neonatal Ventilator

www.wilamed.com

Non-invasive treatment for
the most delicate patients.



CE 0197

Infant nCPAP Ventilation redefined

A new generation in Infant nCPAP Ventilation

WILAflow Elite is a microprocessor controlled, non-invasive infant ventilator, providing most advanced and diversified non-invasive nasal ventilation modes, including apnea wakeup function and automatic leakage compensation.

Direct pressure setting

WILAflow Elite directly sets the value of pressure for fully automatic pressure control.

Direct Oxygen concentration setting

WILAflow Elite uses electronic air/oxygen blender technology with one simple button to set precisely the value of the needed oxygen concentration and can auto-proportionate oxygen and airflow. High-precision flow sensor and proportional valve equipped enable real-time feedback and oxygen concentration precision with +3%.



Precise levels of Oxygen

WILAflow Elite delivers precise levels of O₂ at positive pressure, which helps to keep alveoli open and thus improving oxygenation, while the infant breathes spontaneously. The variable WILAflow CPAP Generator helps to reduce the imposed WOB (work of breathing) during inhalation and exhalation.

Accurate and safe fresh gas delivery by iFlow

The Intelligent Closed-Loop Control System (iFlow) was designed to protect the most fragile patients. iFlow intelligently adjusts fresh gas flow and airway pressure in a closed loop. Proximal pressure monitoring (under the nose) and real-time leakage compensation enable stable pressure output. In case of leakage, iFlow will compensate gas in real time to guarantee stable positive airway pressure. It can compensate leakage up to a maximum of 25%, which is incomparable by traditional CPAP devices.

Safe ventilation weaning

WILAflow Elite features SNIPPV / NIPPV, NCPAP and HFNC modes for safe ventilation weaning.



WILAflow Elite

WILAflow Elite – especially designed for the non-invasive ventilation of preterm infants or newborn infants predisposed with lung disease.

Safe and reliable

WILAflow Elite allows proximal pressure monitoring without being affected by mechanical dead space in the closed loop and compliance. The device accurately measures patient's airway pressure. This is the most recognizable method in the industry.

Optional Abdominal Respiratory Sensor available

The respiratory abdominal sensor enables the clinician to monitor for apnea/low breath rate in both nCPAP and BiPhasic modes. The accessories include the reusable transducer and single-patient-use abdominal sensor. In the BiPhasic trigger mode, the respiratory abdominal sensor and transducer allow patient-triggered pressure assists with breath rate monitoring.



“Adequate humidification is essential to maintain airway clearance, optimize ventilation and improve patient comfort.”

Heated humidification is recommended for nCPAP therapy

The normal functions of the nose and air passages of the respiratory tract are to warm, moisten and filter the inhaled gases before they reach the lungs. In normal respiration, the nasal mucosa and upper airways provide 75% of the heat and moisture supplied to the smaller airways and alveoli. By the time air reaches the alveoli, the inspired gas warms to 37°C at 100% relative humidity (RH).

With nCPAP, the upper airways are not bypassed, but the high gas flows may be drying to the airways, especially to a neonate's underdeveloped lung.

WILAflow Elite

Parameter	Adjustable Range	Step	Ventilation Mode
CPAP	1–13 cmH ₂ O	1–3 cmH ₂ O; 0,2 cmH ₂ O 3–13 cmH ₂ O; 0,5 cmH ₂ O	nCPAP
PEEP	1–13 cmH ₂ O	1–3 cmH ₂ O; 0,2 cmH ₂ O 3–13 cmH ₂ O; 0,5 cmH ₂ O	NIPPV SNIPPV
Pinsp	3–20 cmH ₂ O	3–8 cmH ₂ O; 0,5 cmH ₂ O 8–15 cmH ₂ O; 1 cmH ₂ O	NIPPV SNIPPV
Papnea	3–20 cmH ₂ O	3–8 cmH ₂ O; 0,5 cmH ₂ O 8–15 cmH ₂ O; 1 cmH ₂ O	nCPAP
Manual Ventilation	3–20 cmH ₂ O	3–8 cmH ₂ O; 0,5 cmH ₂ O 8–20 cmH ₂ O; 1 cmH ₂ O	NIPPV SNIPPV nCPAP
	3–25 L/min.	3–10 L/min; 0,5 L/min 10–25 L/min; 1 L/min	HFNC
Manual ventilation time	1–20 s	1 s	nCPAP NIPPV SNIPPV HFNC
Flow	0,5–20 L/min.	0,5–2 L/min; 0,1 L/min 2–10 L/min; 0,5 L/min 10–20 L/min; 1 L/min	HFNC
O ₂ %	21–100 %	1 %	NCPAP NIPPV
Flush O ₂	23–100 %	1 %	NCPAP NIPPV SNIPPV HFNC
Flush oxygen ventilation time	30 - 120 s	30 s	nCPAP NIPPV SNIPPV HFNC
Ti	0,1–20 s.	0,01 s	nCPAP NIPPV SNIPPV HFNC

Parameters and Ranges

Parameter	Adjustable Range	Step	Ventilation Mode
Respiratory Rate	1–120 bpm	1 bpm	NIPPV
Tapnea	OFF; 10–30 s.	5 s	nCPAP SNIPPV
Rb	1–120 bpm	1 bpm	SNIPPV
Pressure		280–600 kPa	
Fresh gas Flow		40 L/min	
Monitoring			
O ₂ %		0–100 %	
CPAP/EPAP Ppeak Pmean		-10–100 cmH ₂ O	
I:E		1:1–1:10	
Te		0,4–30 s	
Rsp		0–200 bpm	
Flow		0–40 L/min	



WILAflow Elite

Part Number	Item	Box / Qty.
101300	WILAflow Elite Infant CPAP Ventilator; including Vent Cart	1
100663	O ₂ low pressure connecting tube, angle plug DIN, 3m	1
100664	O ₂ low pressure connecting tube, angle plug DIN, 5m	1
100666	AIR low pressure connecting tube, angle plug DIN, 3m	1
100667	AIR low pressure connecting tube, angle plug DIN, 5m	1
101302	Signal Box/IEC incl. 2 Body Sensors	1
101301	Body Sensor VX010 100CM	5
101304	Hose Arm Assembly	1
101554	Supply Rack for Sterile Water (50cm)	1
Respiratory Humidification		
101200	AIRcon Gen2 Humidifier, 230V	1
270777	BTS1204A Single limb system for nCPAP, heated (i), A-adapter, autofill chamber, neo-adapter, nCPAP generator, connecting tube for incubator (120cm, 10mm Ø) Connecting tube 40cm	10
CPAP Accessories		
300710	WILAflow CPAP Mask, Size S	10
300711	WILAflow CPAP Mask, Size M	10
300708	WILAflow CPAP Mask, Size L	10
300709	WILAflow CPAP Mask, Size XL	10
300712	WILAflow Nasal Prong, Size S	10
300713	WILAflow Nasal Prong, Size M	10
300714	WILAflow Nasal Prong, Size L	10

Order Information

Part Number	Item	Box / Qty.
300720	WILAbonnet, White, Size 000, 18–20 cm	10
300721	WILAbonnet, Grey, Size 00, 20–22 cm	10
300722	WILAbonnet, Pink, Size 0, 22–24 cm	10
300723	WILAbonnet, Brown, Size 1, 24–26 cm	10
300724	WILAbonnet, Yellow, Size 2, 26–28 cm	10
300725	WILAbonnet, Blue, Size 3, 28–30 cm	10
300726	WILAbonnet, Lite Orange, Size 4, 30–32 cm	10
300727	WILAbonnet, Green, Size 5, 32–34 cm	10
300728	WILAbonnet, Red, Size 6, 34–36 cm	10
300729	WILAbonnet, Orange, Size 7, 36–38 cm	10
300730	WILAbonnet, Turquoise, Size 8, 38–40 cm	10
300731	WILAbonnet, Dark Blue, Size 9, 40–42 cm	10





WILamed GmbH

Medizinische Geräte und Zubehör

Gewerbepark Barthelmesaurach
Aurachhöhe 5-7
91126 Kammerstein (Germany)

Phone: +49 9178 996999-0

Fax: +49 9178 996778

info@wilamed.com

www.wilamed.com

